







HUMIDITY SEPARATORS CLEAN STEAM CENTRIFUGAL SEPARATOR S11

DESCRIPTION

When wet steam is used in sterilization, moisture in suspension reduces the heat transfer efficiency, and the validity of the sterilization process can be compromised.

S11 series centrifugal separators remove moisture from steam pipelines. Steam passes through the separator and, as a result of centrifugal forces, impact and swirling effects, the particles with a heavier specific gravity are separated, such as water droplets and moisture in suspension.

The condensate collected at the bottom of the separator must be automatically drained by a suitable steam trap.



316L stainless steel construction.

No moving parts.

Self draining design.



Internal wetted parts: \leq 0,51 micron Ra – SF1. External: Satin bead blast finish – 1,6 micron Ra.

Other surface conditions see IS PV20.00 E - Technical information.

OPTIONS: Vent connection.

Different kinds of connections and dimensions.

USE: Steam, compressed air and other gases.

AVAILABLE

MODELS: S11 – inline connections.

S11A - S11 with air vent connection.

SIZES: 1/2" to 2".

CONNECTIONS: ASME BPE clamp ferrules.

Others on request.

PACKAGING: Assembling and packaging in a clean room

certified according to ISO 14644-1.

The product is end capped and sealed with recyclable thermo-shrinkable plastic film, to

avoid contamination.

INSTALLATION: Always with the condensate discharge pointing

downwards.

See IMI - Installation and maintenance

instructions.

HOW TO SELECT: Generally, in an existing plant it is advisable to

fit a separator with the same size of the pipe line. Pressure drop is normally negligible. For approximate pressure drop calculation please

consult.





BODY LIMITING CONDITIONS				
ALLOWABLE PRESSURE	RELATED TEMPERATURE			
10 bar	50 °C			
8 bar *	175 °C			
7,4 bar	200 °C			

* PMO – Max. operating pressure for saturated steam. Minimum operating temperature: -10 °C.


Design code: AD-Merkblatt.

CE MARKING – GROUP 2 (PED – European Directive)				
PN 10	Category			
1/2" to 1"	SEP			
11/2" to 2"	1			







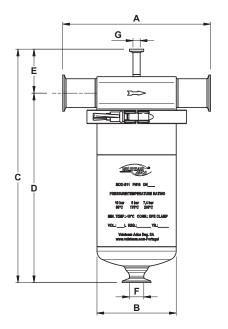



	DIMENSIONS (mm)												
SIZE	Α	В	С	D	E	F1	F2	F3	H1	H2	Н3	VOL. (dm³)	WEIGHT (kg)
1/2"	215	114	345	283,5	62,5	25	25	50,5	9,4	9,4	22,1	2,84	3,8
3/4"	215	114	345	283,5	62,5	25	25	50,5	15,75	9,4	22,1	2,87	3,9
1"	215	114	345	283,5	62,5	50,5	25	50,5	22,1	9,4	22,1	2,9	4,2
11/2"	235	141	416	338,5	77,5	50,5	25	50,5	34,8	9,4	22,1	5,82	7,25
2"	260	141	416	338,5	77,5	64	25	50,5	47,5	9,4	22,1	5,93	7,28

MATERIALS					
POS. Nº	DESIGNATION	MATERIAL			
1	Body	AISI 316L / 1.4404			
2	Cover	AISI 316L / 1.4404			
3	* Clamp	AISI 316L / 1.4404			
4	* Seal	FKM / PTFE			
5	Internals	AISI 316L / 1.4404			

^{*} Available spare parts.

EN 10204 3.1 certificate available if requested along with the order. Remarks: FDA / USP Class VI seals certificate on request. All separators have a serial number. In case of non-standard separators, this number must be supplied if spare parts are ordered.

